

GFS Ultimate is a suspensory fixation device recommended for use in the fixation of ligament or tendon repairs. It combines the convenience of adjustability with the strength and stiffness of conventional fixed-length loop devices.

Features & Benefits

Adjustable loop (Max:60mm to Min:12mm)

- Simplifies surgical procedure
- Eliminates the need to stock multiple sizes
- Can be both tightened and loosened

Loop made of Ultra High Molecular Weight Polyethylene (UHMWPE)

- High strength
- Non-absorbable

Button made of 6Al-4V ELI titanium alloy


- Provides the strength and stiffness required for successful ligament reconstruction

Designed to be used with various surgical techniques/ instrumentation

- Ideal for both transtibial and anterior medial approaches
- Ideal for both outside-in instrumentation and reverse cutting tunnel instrumentation
- Ideal for stepped femoral tunnels or a single diameter femoral tunnel

GFS Ultimate					
Part #	Width of Button (mm)	Length of Button (mm)	Description	Length of Loop (mm)	QTY/Box
10972	3.9	12	GFS, ultimate mini	Max: 60, Min: 12	1
10974	3.9	16.5	GFS, ultimate large	Max: 60, Min: 12	1

Instrumentation for GFS Ultimate					
Part #	Diameter (mm)	Joint	Description	Sterile/ Unsterile	Single-use/ Reusable
10057	-	Knee	GFS Depth Gauge	Unsterile	Reusable
10613U	5	Knee	Headed reamer, w/2.5mm cannulation ¹	Unsterile	Reusable

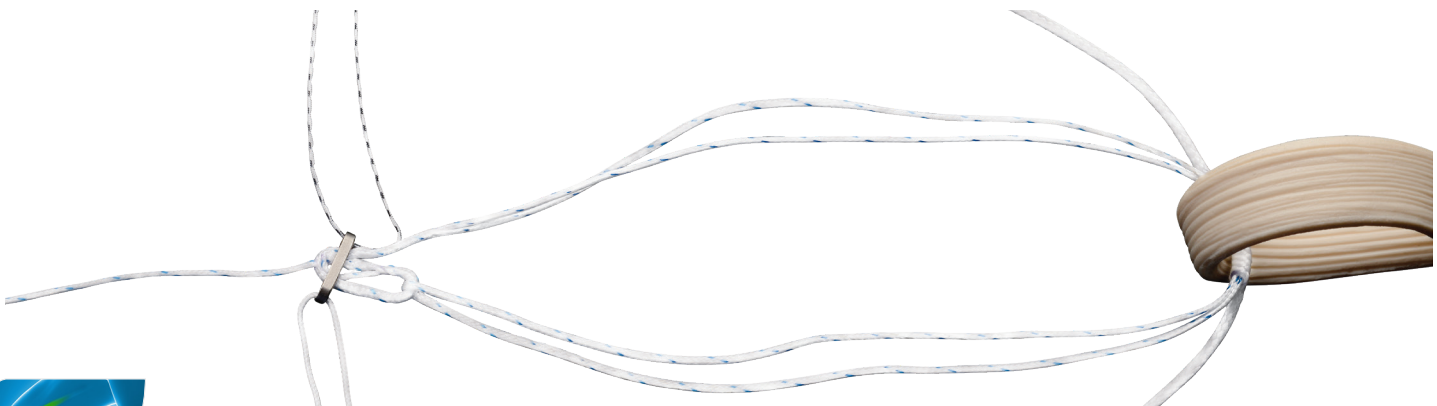

Knee

¹Use w/GFS mini/ultimate